



A100-4A

ELECTRIC OPERATING TABLE

Registration code: SKL-A



SAIKANG Medical

Healthier
your life

A100-4A

FEATURES AND HIGHLIGHTS

Technical parameters

▪ Length/Width:	2115 × 550 mm ± 10 mm
▪ Height Elevating/Lowering:	710-1010 mm ± 10 mm
▪ Forward Tilt/Backward Tilt:	24° ± 3° / 24° ± 3°
▪ Left Tilt/Right Tilt:	23° ± 3° / 23° ± 3°
▪ Backrest (Up/Down):	80° ± 3° / 20° ± 3°
▪ Head Section (Up/Down):	55° ± 3° / 90° ± 3°
▪ Leg Section (Up/Down):	15° ± 3° / 90° ± 3°
▪ Leg Section Outward:	90° ± 3°
▪ Kidney Bridge:	140 mm ± 5 mm
▪ Horizontal Sliding:	350 mm ± 5 mm

Technical configuration

☑ Anesthesia Screen	1 unit
☑ Shoulder Supports	1 pair
☑ Arm Rests (with straps)	1 pair
☑ Lumbar Supports	1 pair
☑ Leg Supports (with straps)	1 pair
☑ Foot Plates	1 pair



Mattress

High frequency heat sealing technology.
Large area of velcro coverage provides strong adhesion.

Hand Controller Storage

All Function Controlled By Hand Controller and Controller in the Column, And Can Be Locked.

Battery

Battery could support about 80-120 complete operation cycles.
Working time could be last for one week according to useage and help to keep the operating table running in the event of power failure.

Kidney Bridge

Manual kidney bridge, operate by head side rocker, more convenience.

Anesthesia Screen Frame

Equipped with a stainless steel anesthesia screen frame for hanging anesthesia curtains, used to shield and cover the patient's head during use.

Shoulder Support

Polyurethane support frame, adjustable in angle, used to stabilize the patient's body.

Body Support

Polyurethane support frame, adjustable in angle, used to stabilize the patient's body.

Arm Rest

Polyurethane widened armrest board with attached straps, providing stable performance and versatile applications.

Leg Rest

Manual multi-angle adjustment for easy use.

Foot Plate

Be used to support patient's feet.

Professional functions

Lowest position of table top

710_{mm}

Highest position of table top

1010_{mm}

Back plate up

0°-80°

Back plate down

0°-20°

Reversed trendelenburg

0°-24°

Trendelenburg

0°-24°

Left lateral tilt

0°-23°

Right lateral tilt

0°-23°

Leg plate up

0°-15°

Leg plate down

0°-90°

Leg plate outwards

0°-90°

Kidney bridge elevation

0-140_{mm}

Head plate up

0°-55°

Head plate down


0°-90°

Horizontal sliding

0-350_{mm}



CONSTANT INNOVATION

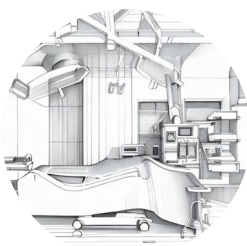


Design



Research



Manufacturing

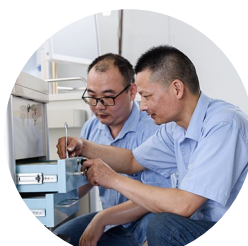

Maintenance



Service

WORLDWIDE BRANCHES

European Branch

Saikang Medical Portugal Office

Av. da Liberdade, nº 101
8150-101
São Brás de Alportel
Portugal
TEL: +351 289 090 212

Russian Branch

LLC 'SAIKANG'

454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

South America Branch

Saikang Medical Technology, INC

Rua Lauro Muller 116, Sala
1503 - Botafogo,
Rio de Janeiro -
RJ, 22290-160
E-mail: mehecobrasil@
gmail.com

East Africa Branch

Saikang Medical Technology Ltd

Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

SAIKANG Medical (Thailand) Co., Ltd.

Ample Tower 909 Room 11/3
Debaratana Road,
Bangna Nuea, Bangkok
10260 Thailand
TEL: +66-936-569469

US Branch

Healthward Technology, INC

17800 castleton st ste 180,
Los Angeles
TEL: +1 (603) 416-3063

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue