



A302

ELECTRIC OPERATING TABLE

Registration code: SKL-A



SAIKANG / Medical

Healthier
your life

A302

FEATURES AND HIGHLIGHTS



Technical parameters

▪ Length/Width:	2000 × 555 mm ± 10 mm
▪ Height Elevating/Lowering:	690-990 mm ± 10 mm
▪ Forward Tilt/Backward Tilt:	24° ± 3° / 24° ± 3°
▪ Left Tilt/Right Tilt:	23° ± 3° / 23° ± 3°
▪ Backrest (Up/Down):	68° ± 3° / 22° ± 3°
▪ Head Section (Up/Down):	33° ± 3° / 90° ± 3°
▪ Leg Section (Up/Down):	18° ± 3° / 85° ± 3°
▪ Leg Section Outward:	90° ± 3°
▪ Kidney Bridge:	110 mm ± 5 mm
▪ Horizontal Sliding:	350 mm ± 5 mm
▪ One-Button Flexion:	220° ± 3°
▪ Reverse Flexion:	100° ± 3°

Technical configuration

<input checked="" type="checkbox"/> Anesthesia Screen:	1 unit
<input checked="" type="checkbox"/> Shoulder Supports:	1 pair
<input checked="" type="checkbox"/> Arm Rests:	1 pair
<input checked="" type="checkbox"/> Lumbar Supports:	1 pair
<input checked="" type="checkbox"/> Leg Supports:	1 pair
<input checked="" type="checkbox"/> Foot Plates:	1 pair

Material

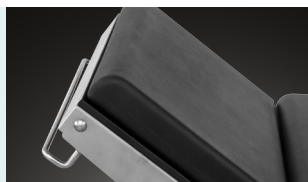
The base material for the table is high thickness 304 stainless steel, seamless welding by robot, with modern appearance and more stable.

**Mattress**

Polyurethane foam mattress, waterproof and seamless, avoid liquid in. The connection between table and mattress is buckle joint, easy to disassemble and install, easy to clean.

**Head plate**

Head plate is with imported gas spring for easy control.

**Kidney bridge**

Manual kidney bridge, operate by head side rocker, more convenience.

**Horizontal sliding**

Extra horizontal sliding, distance 350mm high strengthen phenolic resin table surface, C-arm and X-ray compatible.

**Dual control system**

All functions controlled by both hand controller or control panel on table base column. With action LCD display on hand controller, easy to operate. With control panel on table base column, in case of emergency hand controller failure.

**Leg plate**

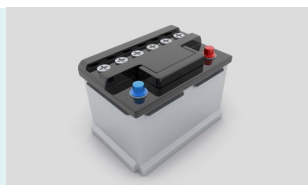
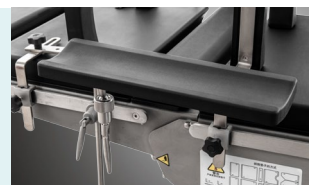
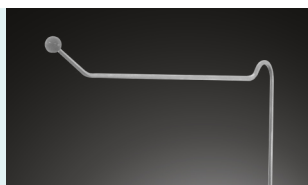
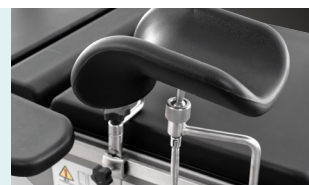
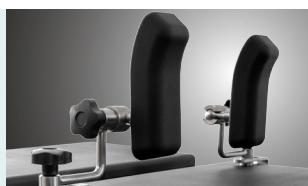
The up and down folding actions of the leg plate are all completed by electric. The leg plate is removable, effortless and smooth.

**Castors**

HAION brand castors, each load bearing is 250kg, smooth and noiseless movement, with longer service life. Pedal brake, easy to move and lock. Unique switch design.

**Battery could support**

Battery could support about 80-120 complete operation cycles. Working time could be last for one week according to useage and help to keep the operating table runing in the event of power failure.

**Armrest****Anesthesia Frame****Leg Crutches****Shoulder Support****Body Supporter**

Professional functions

Lowest position of table top

690mm



Highest position of table top

990mm



Back plate up

0°-68°



Back plate down

0°-22°



Reversed trendelenburg

0°-24°



Trendelenburg

0°-24°



Right lateral tilt

0°-23°



Left lateral tilt

0°-23°



Head plate upward

0°-33°



Head plate downward

0°-90°



Leg section evagination

0°-90°



Leg plate up/down

-18°-85°



Kidney bridge elevation

0-110mm



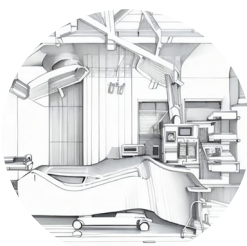
Horizontal sliding

0-350mm





CONSTANT INNOVATION



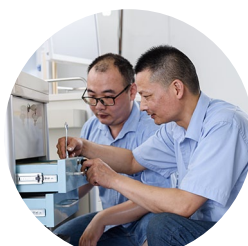
Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

South America Branch

Saikang Medical Technology, Inc
W3400 Corrientes, Argentina
TEL: +55 (11) 2758-0438
Rua Universal, Hollywood,
São Bernardo do Campo,
Brazil
TEL: +55 (11) 2758-0438

East Africa Branch

Saikang Medical Technology Ltd
Mezzanine Floor, Taj Tower,
upper hill, nairobi, Kenya
TEL: +254714084806

South East Asia Branch

Saikang Medical (Philippines) Company Limited
Manila, The Philippines
TEL: +63 953 044 7864

US Branch

Healthward Technology, INC
17800 CASTLETON ST STE 180,
CITY OF INDUSTRY CA 91748
TEL: +1 (603)416-3063

Russia Branch

LLC 'SAIKANG'
454004, Russian Federation,
Chelyabinsk,
Universitetskaya
Naberezhnaya, Russia
TEL: +7 (919)300-00-06

Mid East Branch

Saikang Medical Supplies (Mid East) LLC
Musaffah
P.O Box 53200, Abu Dhabi,
UAE
Tel: +971 585710652

SAIKANG Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China
☎ +86 18021231901 · 🕒 24/7
✉ export@saikangmedical.com
🌐 www.saikangmedical.com



Company video



General catalogue